


ORDERING INFORMATION

ELEGANCE Silicone

is available in the following sizes:

Size	Approx Length	Application	Left	Right	Color Ranges
STANDARD COLORS					
8 ¹ / ₄	645mm	Male	SG300	SG301	E1 to E18
7 ³ / ₄	645mm	Male	SG302	SG303	E1 to E18
7 ¹ / ₄	645mm	Female	SG304	SG305	E1 to E18
7 ¹ / ₄	645mm	Youth	SG306	SG307	E1 to E18
TRUEFINISH™ COLORS					
8 ¹ / ₄	645mm	Male	TSG300	TSG301	E1 to E18
7 ³ / ₄	645mm	Male	TSG302	TSG303	E1 to E18
7 ¹ / ₄	645mm	Female	TSG304	TSG305	E1 to E18
7 ¹ / ₄	645mm	Youth	TSG306	TSG307	E1 to E18

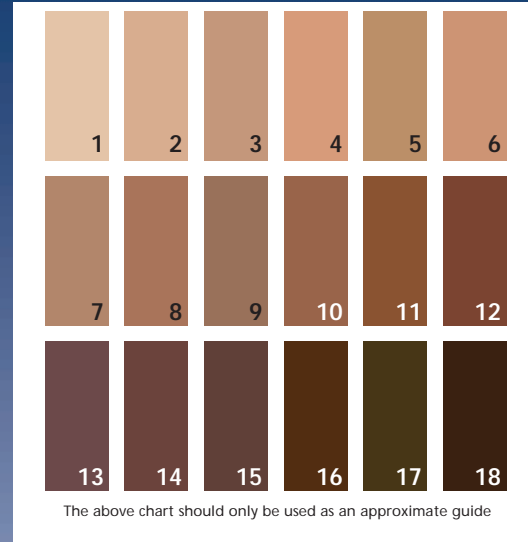

ELEGANCE PVC

is available in the following sizes:

Size	Approx Length	Application	Left	Right	Color Ranges
STANDARD COLORS					
8 ¹ / ₄	645mm	Male	CG300	CG301	E1 to E18
7 ³ / ₄	645mm	Male	CG302	CG303	E1 to E18
7 ¹ / ₄	645mm	Female	CG304	CG305	E1 to E18
7 ¹ / ₄	645mm	Youth	CG306	CG307	E1 to E18
TRUEFINISH™ COLORS					
8 ¹ / ₄	645mm	Male	TG300	TG301	E1 to E18
7 ³ / ₄	645mm	Male	TG302	TG303	E1 to E18
7 ¹ / ₄	645mm	Female	TG304	TG305	E1 to E18
7 ¹ / ₄	645mm	Youth	TG306	TG307	E1 to E18



Available in
18 color shades
and 4 sizes

ORDERING EXAMPLES

TSG 300/5

Glove number followed by color

TG300/5

Glove number followed by color

DISTRIBUTED IN THE USA BY:



Liberating Technologies, Inc.
325 Hopping Brook Road, Suite A
Holliston, MA 01746-1456
Phone. 508 893 6363 Fax. 508 893 9966
Toll free. 800 437 0024
E-mail. info@liberatingtech.com
www.liberatingtech.com

L25006 ISSUE 1
© RSLSTEPPER SEPTEMBER 04



ELEGANCE

Extra Length Silicone & PVC
Cosmetic Gloves

Prosthetic Product Solutions

The first full length range of Silicone & PVC Cosmetic Gloves

ELEGANCE Gloves

Responding to clinical requests, we have made our ELEGANCE range of Silicone and PVC Gloves in our new extra long style. This allows smooth transitions over the proximal socket brim and, where needed an extension to mid upper arm, providing a unique continuous cosmesis.

- Extra length as standard in 4 sizes,
- Suitable for all Adult Otto Bock™ system hands
- Advanced micro-pigmentation finish available
- Suitable for Myo electric and Cable powered hands
- Available in 18 color shades and in two ranges, Standard and Truefinish™

ELEGANCE Silicone

- Tufsil 3™ New silicone with proven toughness
- Easy to clean with soap and water
- Extra stretch minimises the energy used by the hand

Using Tufsil 3™, our latest innovation in polymer technologies, these gloves offer the toughness associated with PVC but provide the ease of cleaning provided by silicones. In addition, Tufsil 3™ allows additional stretch, reducing the energy expended by the hand motor each time it opens. Properties of this polymer include improved tear resistance proven by clinical evaluations and up to 15% stretch, to allow even bulbous sockets to be accommodated. These gloves are easy to fit and remove, allowing quick access to the hardware for maintenance and adjustment. They are simple to keep clean using soap and water.

ELEGANCE gloves are available in a standard finish, as well as Truefinish™ (incorporating our micro pigmentation feature for added realism and depth). Both finishes are supplied with hand painted nails.




ELEGANCE PVC

These gloves provide a cost-effective solution that will cover prostheses fitted with adult Otto Bock™ electric and system hands. The additional length of these gloves allows them to extend up to and over the proximal edge of the socket, providing a smooth transition to the skin. For trans-humeral and elbow disarticulation prostheses the glove can be trimmed at the elbow fold and the remaining proximal portion used to form a discontinuous cosmesis.



Details of palm and skin texture

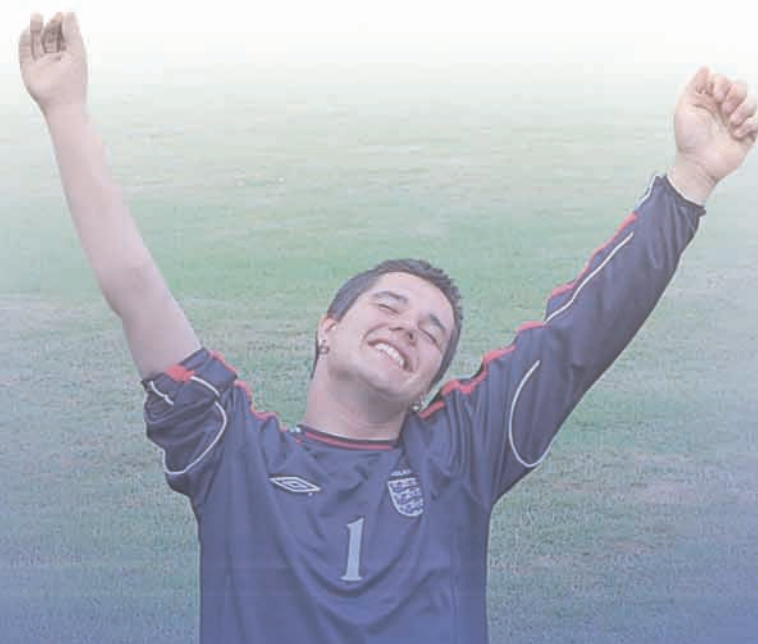

ELEGANCE
ELEGANCE
ELEGANCE
ELEGANCE